

OTOCLEAR® EAR IRRIGATION SYSTEM

Kit includes:
1 spray wash bottle
20 otoclear tips
1 ear basin



25858

- 25858 OTOCLEAR EAR SPRAY WASH
- manual, intermittent lavage
- delivers greater volume compared to syringes or other manual devices
- efficient design allows for two handed procedure




Thermometer for water temperature



Improved patient comfort thanks to the soft "gentle touch" design

- 25861 OTOCLEAR EAR IRRIGATION TIPS
single use - box of 40 - spare
Flared design protects against over-insertion. The tip directs three streams of water towards the walls of the ear canal, away from the tympanic membrane, creating a turbulent lavage that effectively dislodges the cerumen. The exit portal controls backsplash and directs drainage into a basin.
- 25862 EAR IRRIGATION BASINS - box of 3 - spare

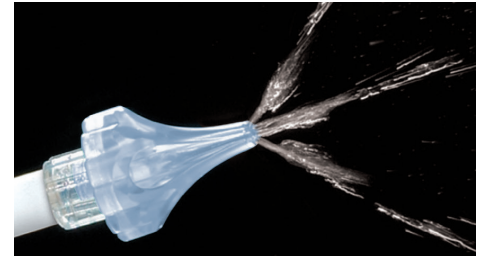

Patented 3-stream action

LOOP-SPOON EAR CURETTES

- 31470 GIMA EAR CURETTES - box of 10 pcs
Double function currettes: one side spoon, one side loop.
Each curette is individually packed in a single pouch.

31470



Loop side

Spoon side

BIONIX EAR CURETTES



FLEXLOOP:
For normal curettage
Most flexible of all
Least traumatic
Rounded tips for comfort
4 mm wide tip



MICROLOOP:
For compact cerumen
More rigid construction
Tip can be angled for added strength
Rounded tip for comfort
3 mm wide tip



INFANTSCOOP:
For the smallest ears
For when a small scoop works better
Can be used through an otoscope • Flexible
2 mm wide tip



CERASPOON:
For viscous cerumen
For extensive cerumen buildup
For when a large scoop works better
4 mm wide tip

The complete cerumen management system

- 31462 FLEXLOOP - white - box of 50 pcs
- 31464 INFANTSCOOP - light blue - box of 50 pcs
- 31460 MICROLOOP - green - box of 50 pcs
- 31467 CERASPOON - yellow - box of 50 pcs




Our four styles of currettes give you the most effective, economical and safe way to remove cerumen - whatever the size of ear or consistency of cerumen.

Patient comfort: Metal currettes are cold, hard and sharp. BIONIX ear currettes with rounded plastic tips reduce trauma and risk of laceration. Latex free.

Effectiveness: Ear currettes can duplicate the effectiveness of metal currettes in almost any situation and their polypropylene composition makes them virtually unbreakable.

Disposability: Ear currettes don't require expensive and time-consuming sterilization and eliminate the risk of cross-contamination.

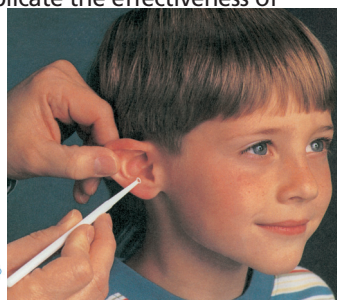

LIGHTED EAR CURETTES

- 31450 KIT MIXED LIGHTED CURETTES
Kit is composed by 50 lighted currettes (8 FlexLoop and 7 of each of below other 6 tip styles) + LED illuminator + magnification lens



MicroLoop®
3 mm diameter loop



FlexLoop®
4 mm diameter tip



InfantScoop®
2 mm diameter tip



AngleLoop®
Angled 4 mm diameter loop



WaveCurette™
new wave like 2.5 mm tip



VersaLoop®
3 mm wide tip



CeraSpoon®
4 mm wide tip



Scan QR code with your web enabled camera phone to see Video



31458

Lighted forceps for foreign body removal with LED illuminator



31459

Lighted ear curette 3 mm tip with LED illuminator

- 31461 KIT MICROLOOP
- 31463 KIT FLEXLOOP
- 31465 KIT INFANTSCOOP
- 31466 KIT ANGLELOOP
- 31468 KIT VERSALOOP

All kits are composed by 50 lighted currettes + LED light source + magnification lens

- 31458 LIGHTED FORCEPS - box of 10 + magnification lens + LED illuminator
- 31459 LIGHTED ARTICULATING CURETTES - box of 25 + magnification lens + LED light source

Ear curette tip has a 90° range of motion. It articulates to assist in removal of cerumen from the ear canal

