

# VETERINARY

## ENDOSCOPE SERIES



**HUGER**  
ENDOSCOPY SYSTEM

# AGVE 68 HD SERIES

- High definition system, CMOS technology with resolution as 1280\*720p, frame rate as 60fp/s and IPS display technology
- High quality still and dynamic images
- New endoscope handle of ergonomics design
- Compact design, integrated with processor, light source and monitor in one unit



## AGVE-68 Series Endoscope



Type	AGVE-68HQ	AGVE-68HS	AGVE-68HB	AGVE-68HAL	AGVE-68HBL	AVB-68Q/QL	AVB-68PL
Field of View	145°					120°	
Depth of View	3~100mm					3~50mm	
Tip Deflection	Up/Down:210°/90° Left/Right:100°/100°			Up/Down:180°/180° Left/Right:160°/160°		UP180° DOWN130°	UP130° DOWN130°
Distal Rigid Diameter	9.2mm	8.0mm	9.6mm	12.9mm	18.0mm	5.9mm	4.3mm
Insertion Tube Diameter	9.2mm	8.0mm	9.6mm	12.9mm	18.0mm	5.8mm	3.8mm
Working Channel	2.8mm	2.2mm	2.8mm	3.7mm	6.0mm	2.0mm	N.A.
Working Length	1500mm	1500mm	1500mm	3000mm	2500mm	600/1000mm	600mm

Remarks: Special specifications are customizable.

# VIS-68 HD VIDEO CENTER

Compact design with processor, LED light source and 12 inch IPS monitor in one unit, portable and space-saving



VIS-68-001



VIS-68-002

## VIS-68 Video Center

Type	VIS-68
Lamp	LED 40W
Life of Lamp	35000 hours
Color Temperature	5300K
Pump	Durable and quiet pump
Function	Automatic white balance, Color adjustment, Color saturation adjustment, Brightness adjustment, File establishment, Image capturing, Video recording, SD card compatibility
Video Out	DVI*1, HDMI*1
Monitor	12 inch IPS monitor (001), 26 inch display screen (002)



# EC-68 CAMERA HEAD

- High quality CMOS sensor
- Light-weight body, making the scope more maneuverable, reducing stress and fatigue even during long procedures



## EC-68 Camera Head

<b>Size</b>	Camera Head	Φ32mm*70mm
	Cable	Φ5mm*3M
<b>Observation</b>	Pickup System	CMOS sensor of 4 million pixels
<b>Classification (electro-magnetic equipment)</b>	Degree of protection against electric shock	Type BF
	Degree of protection against explosion	This product should be kept away from flammable gasses
<b>Cleaning/Disinfection Sterilization</b>	Cleaning/Disinfection	Wipe with gauze wetted with disinfection solution
	Sterilization	Ethylene oxide gas

Remarks: Adapters are optional.

# VETERINARY RIGID ENDOSCOPE

- Various rigid endoscopes



Urethra-cystoscope, arthroscope, hysteroscope, otoscope (with sheath)



Urethra-cystoscope, arthroscope, hysteroscope, otoscope (without sheath)



Laparoscope

# USB DIGITAL SERIES



- High quality image
- Easy system connection, USB connector to PC
- Easy for carrying, compact endoscope design integrated with LED flashlight



## USB Digital Endoscope

Type	UCB-66A	UCB-66L	UCN-66A
Field of View		110°	
Depth of View		3-100mm	
Tip Deflection	Up/Down:180°/130°		Up/Down:130°/130°
Distal Rigid Diameter		5.9mm	3.8mm
Insertion Tube Diameter		5.9mm	3.8mm
Working Channel		2.0mm	NO
Working Length	600mm	1000mm	300mm

# SOFTWARE IMAGE PROCESSING SOFTWARE

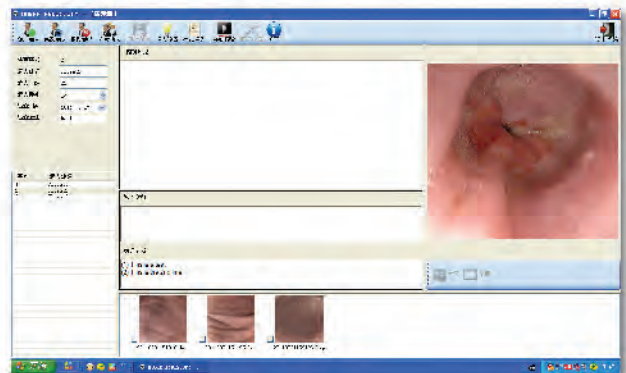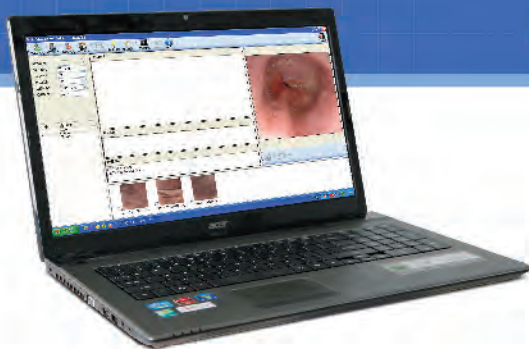

- User friendly and easy to learn
- Automatic white balance
- Color temperature adjustment
- Automatic electronic shutter
- Various image display shape and size options
- High quality images and videos capturing, recording, storage and replay
- Standard report printing
- High compatibility, compatible with Huger video GI endoscope, bronchoscope and ENT scope



# ACCESSORY

- Compatible with GI endoscope and bronchoscope

## Accessory Specification

01



Rat-Tooth Type Grasping Forceps

Alligator Jaws Type Grasping Forceps

02



03



Pelican Jaws Type Grasping Forceps

Rat Tooth with Alligator Jaws Type Grasping Forceps

04



05



Fold Angular Basket Type Grasping Forceps

3-Prong Type Grasping Forceps

06



# ENDOSCOPY TROLLEY

- 4 secure shelves, easy to access to all the equipment in the endoscopy unit
- 1 Adjustable arm, supporting 2 endoscopes
- 1 drawer, convenient for stuff storage
- 4 Castors including 2 lockable ones, easy for trolley moving and locking
- 4 power sockets with central power button and bright LED power-on indicator



**HUGER**  
ENDOSCOPY SYSTEM

**HUGER MEDICAL INSTRUMENT CO., LTD**

Building 26, Xinzhuan Highway 3825, Songjiang District, 201619, Shanghai, China

Tel: +86 21 6762 6235 / 6762 6236

Fax: +86 21 6769 1721

[www.huger.cn](http://www.huger.cn)